


# Flow Rate Sensor Specifications

## Media Isolated Microfluidic Flow Rate Sensor

- Measures up to 80  $\mu\text{L}/\text{min}$
- Non-invasive measurement
- Simple UART Output
- 15 Seconds Response Time

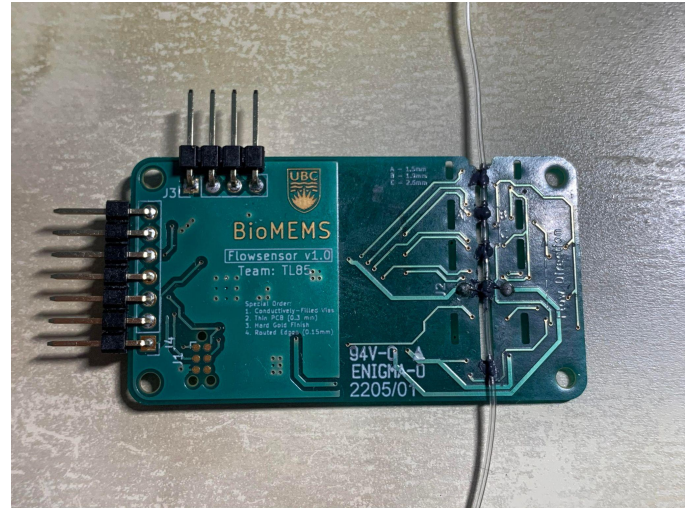

## Electrical Characteristics

Parameter	Min	Nominal	Max	Unit
Supply Voltage	1.8	3.3	3.6	V
Supply Current	-	10	600 (1)	mA
Response Time	14	15	16	s

Note 1: Drains up to 600 mA while applying power to the heater wires

## Measurement Characteristics

Parameter	Min	Nominal	Max	Unit
H2O Flow Rate Full Scale	-	-	80	uL/min
H2O Accuracy Below 55 uL/min	-	-	5	% full scale
H2O Accuracy Below 80 uL/min	-	-	18	% full scale

All test results taken with the following parameters

- Supply Voltage = 3.3V
- Ambient Temperature = 23 C
- Water = DDW (Distilled, Deionized Water)

## Operation & Storage Characteristics

Parameter	Min	Nominal	Max	Unit
Storage Temperature	-10	23	60 (2)	C
Operating Temperature	10 (3)	23	50 (3)	C
Storage & Operating Humidity	-	-	95	%RH
Maximum Reagent temperature	-	-	30 (4)	C

Note 2: The Tygon tubing may have a lower storage temperature than this figure

Note 3: These figures are best guesses, they have not been tested

Note 4: The maximum temperature the reagent will be heated to during the measurement process

## Mechanical Characteristics

Parameter	Min	Nominal	Max	Unit
Measurement Tube Inner Diameter (5)	0.184	-	0.254	mm
H2O Pressure Drop @ full scale	-	-	10 (6)	mbar
Total Internal Volume	-	-	2	uL
Max Pressure	-	-	2000 (7)	mbar


Note 5: Measurement tube consists of alternating sections of stainless steel 29G tubes and 0.01"ID, 1/32"OD Tygon tubes

Note 6: Theoretical, calculated with 80uL/min as flow rate

Note 7: Experimental result where a force equivalent to 2 bars of pressure were applied to the syringe connected to the flow sensor


## Specification Charts

### %Error of Full Scale

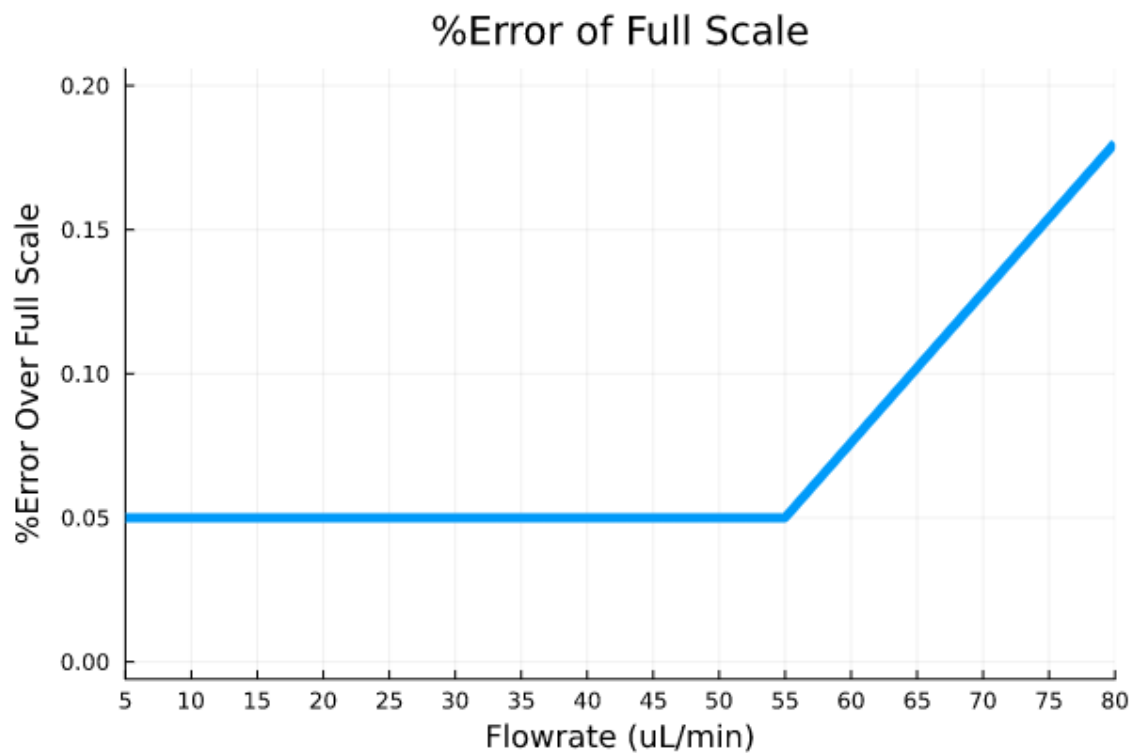


## Benchmark Result with Sensirion SLI-0430

Sensirion SLI-0430 vs Low Cost Flowsensor

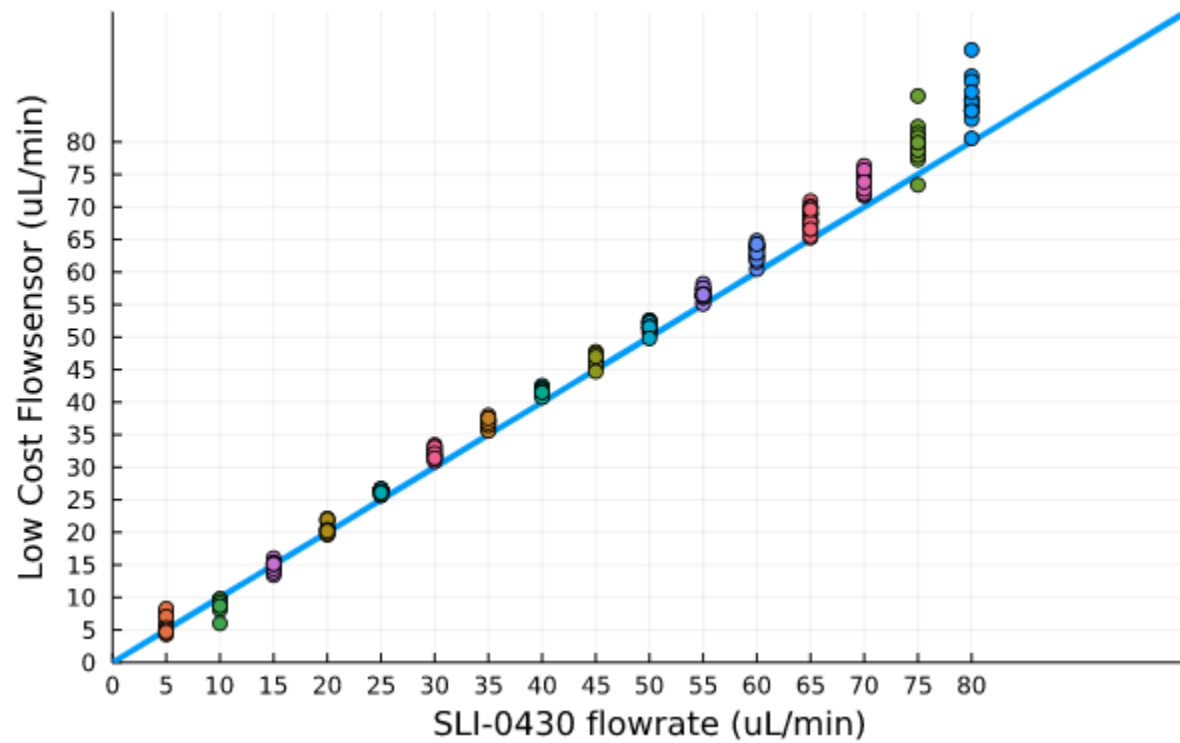

# Pinout

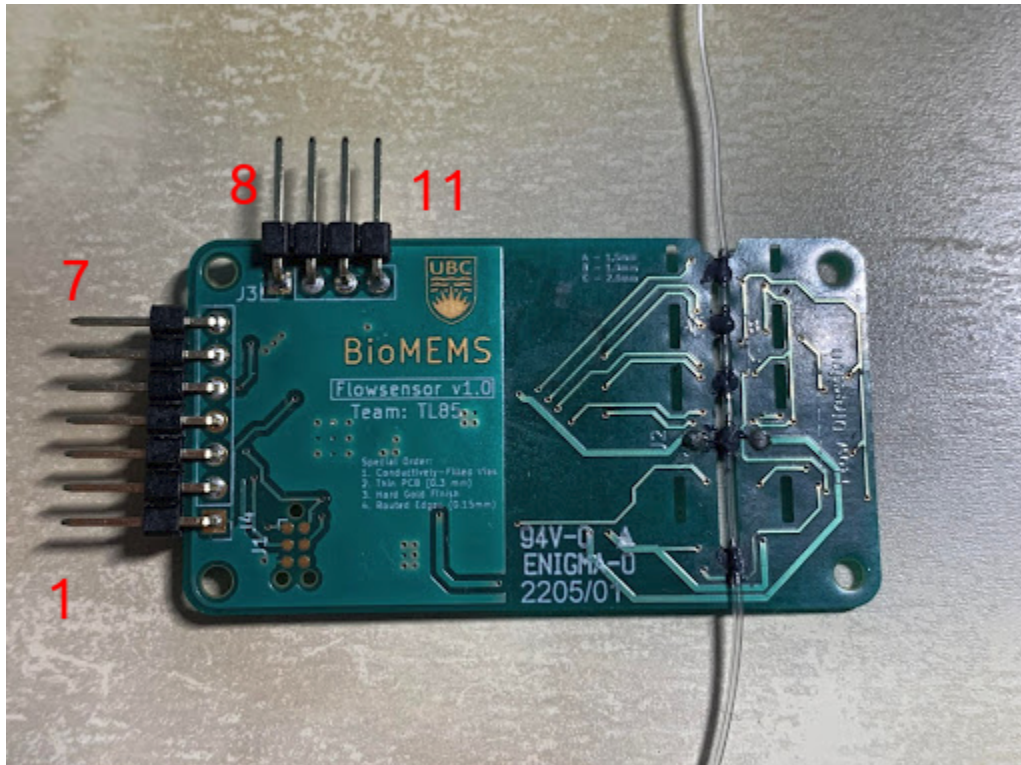

In the picture above, the pins are numbered 1 to 7 from bottom to top and 8 to 11 from left to right.

Pin Number	Pin Description
1, 9	VDD
2	SWCLK
3, 8	GND
4	SWDIO
5	NRST
6	VCP_TX
7	VCP_RX
10,11	Do not connect

# Calibration

Manual calibration is required. Here are the resources required for the calibration step.

1. STLINK-V3
2. Access to the github repo <https://github.com/jackyruth/UBC-BioMEMS-Capstone-Code>
3. A system to drive a reference flow rate

Here are the steps to calibrate the sensor

1. Connect the STLINK-V3 to the sensor
2. Program 'stm32/data\_collection' firmware
3. Drive a constant flow rate through the sensor and store the data
4. Push the data through 'julia/flowrate\_detection.jl' which outputs three calibration constants A, B and C
5. Change the calibration constants in 'julia/functional\_firmware'
6. Program the 'julia/functional\_firmware' and the device should start returning proper flow rate measurements

# Mechanical Specs

Dimensions of the sensor are 28.75x53.00 mm.

The PCB is 0.3mm thin with conductively filled vias to reduce thermal resistance.

Mounting holes are for M2 screws.

The wetted materials are stainless steel tubing cut from 29G needles and 0.01" ID 1/32" OD Tygon tubing.